

Anis Expert Anterior Teeth



Anterior teeth Anis-Expert have exceptional tooth mould and surface texture.

Comprise four layers: Cervical (1), dentin (2), incisal (3) and back incisal (4).

Smooth shade of cervical area gives possibility to leave neck open when needed.

Dentin core of the tooth is encapsulated by incisal on all sides. This gives the tooth a true-to-nature appearance, and makes it possible to grind the lingual area of the tooth without loosing tooth color or getting too translucent incisal edge.



Incisal edge has a smooth transition of layers without having transparent areas. This simplifies the technician's job when setting up and gives highly esthetic results.

Exact color matching to A-D shades. Layering is designed in such a way that each tooth has a "depth" together with shade stability.

Anis Expert Posterior teeth



Anis Expert posterior teeth are easy to set-up while having anatomic-looking mould.


Comprise four layers: Cervical (1), dentin (2) and incisal (3).

Layering is designed in such a way that vestibular and lingual side of the tooth have almost the same shade, which gives a natural-looking effect even when turning the tooth is needed.

Notch on the bottom side reduces grinding in cases with low alveolar height.

Lingual side is high enough to use teeth for combined work.

Larger occlusal area material volume gives more possibilities for grinding.

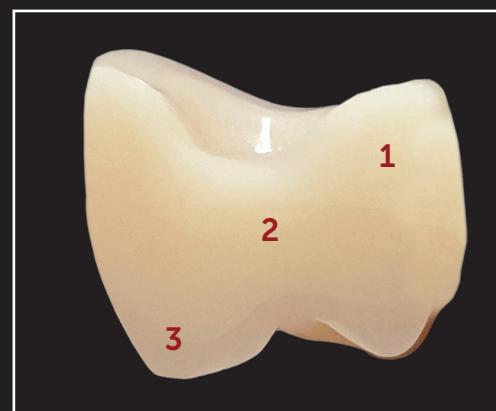

Translucency of occlusal area is reduced to avoid difference in shade and translucency when larger grindings.

High-pressure 3 layer pressing technique of cross-linked VHMV (very high molecular weight) acrylics makes posterior teeth Anis Expert even more stronger.